

High Performance Liquid Chromatography Identification, Quantification, and Fractionation of a Suspect Allergen, 4-Chloro-3-methylphenol, in a LLNA Positive Metalworking Fluid

Kenneth K. Brown, Research Chemist, Division of Applied Research and Technology, National Institute for Occupational Safety and Health, 4676 Columbia Parkway MS-R7 Cincinnati, Ohio 45226 krb1@cdc.gov

Keisha Robinson, NIOSH-IMHOTEP Intern, Howard University, Department of Biology, 415 College Street NW Washington, DC 20059

ABSTRACT

A Metalworking Fluid (MWF) was obtained that produced an allergic reaction in the local lymph node assay (LLNA) with an EC₃ = 4%, the EC₃ being the estimated concentration needed to provoke a 3-fold allergy response. High performance liquid chromatography (HPLC) was used to separate, identify, and isolate the suspected allergen. The biocide, 4-chloro-3-methylphenol, was detected in the MWF as a chromatographic peak matching the retention time of an external standard. The technique of standard addition was used to quantify and confirm the presence of 4-chloro-3-methylphenol at about 1% (w/w). Preparative HPLC was used to fractionate 1 gram of MWF separating the 4-chloro-3-methylphenol fraction from the remaining MWF. The two mobile phase solutions were concentrated back into a MWF and a 4-chloro-3-methylphenol fraction. The original MWF and the reconstituted 4-chloro-3-methylphenol and MWF fractions were also analyzed by gas chromatography - mass spectrometry to confirm the isolation of the biocide.

INTRODUCTION

Estimates indicate that more than 13 million workers in the United States are potentially exposed to chemicals that can be absorbed through the skin. A worker's skin may be exposed to hazardous chemicals through direct contact with contaminated surfaces, deposition of aerosols, immersion, or splashes. When substantial amounts of chemicals are absorbed, systemic toxicity can result. Contact dermatitis can also result when chemicals are absorbed through the worker's skin and is one of the most common chemically induced occupational illnesses, accounting for 10-15 percent of all occupational illnesses at an estimated annual cost of at least \$1 billion. (NIOSH 2006)¹.

Some 1.2 million workers in machine finishing, machine tooling, and other metalworking and metal-forming operations are potentially exposed. Workers can be exposed to the fluids by breathing aerosols generated in the machining process, or through skin contact when they handle parts, tools, and equipment covered with the fluids. The National Institute for Occupational Safety and Health (NIOSH) defines MWF aerosol as the mist and all contaminants in the mist generated during grinding and machining operations involving products from metal and metal substitutes. (NIOSH 2006)¹

Worldwide usage of MWF is estimated to be 2 billion liters per year (Cheng, Phipps, & Alkhaddar, 2005)². MWFs are divided into four classes based on their oil and water content: insoluble, soluble, semi-synthetic, and synthetic. Insoluble (or straight)

oil MWFs function mostly as lubricants, and are not diluted with water. Similar to straight oil MWFs, soluble oil MWFs are used to cool and lubricate tool surfaces, but are diluted with water. Semi-synthetic MWFs contain small amounts of oil and additives, while synthetic MWFs have no oil at all in its formula (NIOSH 1998)³.

Besides water and oil, MWFs can contain hundreds of different compounds and additives such as amines, surfactants, emulsifiers, detergents, and biocides. These intentional additives are chemicals used to modify the formula improving or enhancing its performance (Sheehan, 1999)⁴. Biocides are often times used to destroy the variety of microorganisms or fungi that can survive in MWFs (Cheng et al., 2005²). If gram negative bacteria are present, endotoxins are released. These endotoxins can contaminate MWFs and may enhance the allergic responses in exposed workers (Lim et al.)⁵.

MWF exposure assessment has been a priority research area at the National Institute for Occupational Safety and Health (NIOSH) as evident by research dating back to the early 1970s (Glaser et. al.)⁶. In a 1998 document, Criteria for a Recommended Standard: Occupational Exposure to Metalworking Fluids, NIOSH recommended that the level of exposure for any given worker should not exceed 0.4 mg/m³ of air (thoracic particulate mass) as a time weighted average concentration for up to a ten hour day during a forty hour week. Gravimetric and infrared spectrophotometric techniques, published in the NIOSH Manual of Analytical Methods as Methods 0500, 5024, and 5526, can be used to estimate MWF exposure.⁷ When using gravimetric or

infrared techniques, the results are limited to estimating the mass or concentration level of the MWF without the identification of its chemical composition (Verma et al. and Raynor, Kim, & Bhattacharya)^{8,9}. Because chemicals have very different toxicities, the toxicity of an exposure cannot be assessed without knowledge of the chemical composition.

Little is known about the exact chemical make-up of each MWF because of the competitiveness of the industry and trade secrets. A method that can separate and identify the components of a MWF is needed to assist workers who develop allergic contact dermatitis. Analytical chromatography is a technique used to separate complex mixtures. In HPLC, compounds are separated on a chromatography column and elute off as purified components. HPLC has been utilized to identify the presence of contact allergens and to isolate allergen bands or fractions (Lee, Ong, Koh, & Ong, and Wahl, Meineke, & Maasch).^{10,11}

In this project, semi-preparative HPLC was used to identify and isolate 4-chloro-3-methylphenol from a metalworking fluid. About 1000 mg of sample is needed for a Local Lymph Node Assay (LLNA). 4-chloro-3-methylphenol has also been used in skin cosmetics. Andersen et. al.¹² determined that 4-chloro-3-methylphenol was a sensitizer in guinea pigs using 5 topical preparations and the cumulative contact enhancement test. However, that test may have overestimated the sensitizing potential of 4-chloro-3-methylphenol because later studies using human skin patch testing showed only 2% of the 1462 subjects were allergic to this agent. They concluded that the results from guinea pig allergy tests cannot stand alone, but have to be validated by other sources of

information. Later in 1997, an expert panel reported in the International Journal of Toxicology a “Final Report on the Safety Assessment of p-Chloro-m-Cresol” that concluded the available data was insufficient to support the safety of 4-chloro-3-methylphenol for use in cosmetic products.¹³ As with cosmetics, the safety assessment of 4-chloro-3-methylphenol in metalworking fluids that contact the skin has insufficient data, and thus this study is to provide more data for safety assessment in metalworking fluids.

2 METHODS

2.1 Instrumentation

The HPLC system consisted of HP-1050 modular units that included an injector, a UV detector, and a quaternary solvent pump (Agilent Technologies, Palo Alto, CA, USA). UV light adsorption was monitored at 254 nm for analyte detection. The column was a 10 mm X 300 mm XTerra® Prep MS-C18 (Waters, Milford, MA, USA). The column had a 10 µm particle size, a 0.65 cm³/g pore volume, and an average pore diameter of 113 Å. Mobile phase flow rate was 5.0 mL/min at room temperature, and the injection volume was 100 µL. This method used a ternary mobile phase system, three solutions. The Dionex AI-450 Chromatography Data Acquisition Software was run on a Microsoft Windows PC system (Dionex Corp., Sunnyvale, CA, USA).

2.2 Chemicals

The mobile phase solvents were ACS-HPLC grade hexane, 2-propanol, and acetonitrile from Burdick & Jackson (Honeywell Corp., Morristown, NJ, USA). Mobile Phase A used purified water from a Milli QTM Nanopure Water System Model D4751 (Barnstead Int., Dubuque, Iowa USA). Formic acid, p.a. (pro analysi) grade, was added to mobile phase A at 0.1%, v/v, (Acros Chemicals, New Jersey, USA). Mobile Phase B was 33.3% 2-propanol and 66.6% acetonitrile. Mobile Phase C was 33.3% 2-propanol and 66.6% hexane. A "universal solvent mixture" (USM) solution of methanol, 2-propanol, and hexane (1:1:1 v/v/v) was prepared to dissolve the MWF concentrate.

2.3 Identification and Quantification of 4-chloro-3-methylphenol in MWF

A MWF sample was obtained from the National Toxicology Program, NTP. A standard of 4-chloro-3-methylphenol (chlorocresol) was obtained from Chem Service, Inc. (West Chester, PA, USA). The presence of 4-chloro-3-methylphenol in the MWF was determined by matching chromatographic peak retention times of the samples against the external standards. Then, the standard addition technique was performed to help confirm the identity of the peak and estimate the original amount of 4-chloro-3-methylphenol in the MWF. Four solutions were prepared for standard addition analysis; A, B, C, and D. Solution (A) consisted of only 4-chloro-3-methylphenol external standard, 7.80 mg diluted in 10.00 mL of USM (0.780 mg/mL). Solution (B) consisted of only MWF, 1.878 g diluted in 25.00 mL of USM (75.12 mg/mL). This MWF solution was sonicated at 55°C and its density was determined. In Solution (C), 1.23 mg of 4-chloro-3-methylphenol standard was added to 1.50 mL of MWF solution (+0.82 mg/mL).

Solution (D) consisted of 5.3 mg of 4-chloro-3-methylphenol added to 1.50 mL of MWF solution (+3.53 mg/mL). Using HPLC analysis, the retention times and peak areas of the 4-chloro-3-methylphenol were determined for each solution. Injection volume was 100 μ L(0.1 mL). Table 2.3 was used for the ternary mobile phase gradient program

| Table 2.3 Mobile phase gradient program for analytical liquid chromatography. | | | | |
|---|------------------|------------------|------------------|-------------|
| Time (minutes) | % Mobile Phase A | % Mobile Phase B | % Mobile Phase C | Flow mL/min |
| 0 | 100 | 0 | 0 | 0.1 |
| 0.1 | 100 | 0 | 0 | 5.0 |
| 10.0 | 50 | 50 | 0 | 5.0 |
| 60.0 | 0 | 100 | 0 | 5.0 |
| 90.0 | 0 | 0 | 100 | 5.0 |
| 90.1 | 0 | 100 | 0 | 5.0 |
| 100 | 0 | 100 | 0 | 5.0 |
| 100.1 | 100 | 0 | 0 | 5.0 |
| 110 | 100 | 0 | 0 | 5.0 |
| | | | | |

2.4 Fractionation Procedure

To fractionate the 4-chloro-3-methylphenol from the MWF, the HPLC method was used, but with a truncated mobile phase program to save time as found in Table 2.4. The injector was programmed for multiple 100 μ L injections of MWF solution, ten times for a total of 1.0 mL of solution, onto the column in order to deliver 75 mg of MWF. This was done by programming 10 injections per vial with a run stop time of 0.1 minutes for the first 9 injections, and then manually changing the run stop time to 60 minutes after the ninth injection. The HPLC effluent containing 4-chloro-3-methylphenol was diverted and captured during a predetermined retention time window into a 10 mL test tube. It took 17 HPLC runs to process 1,275 mg of MWF into two fractions, a 4-chloro-3-

methylphenol fraction and a MWF fraction without 4-chloro-3-methylphenol. To determine the efficiency of this method, cool-on-column gas chromatography - mass spectrometry (COC-GC-MS) analysis was done on each fraction for 4-chloro-3-methylphenol.

| Table 2.4 Mobile phase gradient program for semi preparative liquid chromatography. | | | | |
|---|------------------|------------------|------------------|-------------|
| Time (minutes) | % Mobile Phase A | % Mobile Phase B | % Mobile Phase C | Flow mL/min |
| 0 | 100 | 0 | 0 | 0.1 |
| 0.1 | 100 | 0 | 0 | 5.0 |
| 10.0 | 50 | 50 | 0 | 5.0 |
| 20.0 | 40 | 60 | 0 | 5.0 |
| 20.1 | 0 | 100 | 0 | 5.0 |
| 30.0 | 0 | 100 | 0 | 5.0 |
| 30.1 | 0 | 0 | 100 | 5.0 |
| 40.0 | 0 | 0 | 100 | 5.0 |
| 40.1 | 0 | 100 | 0 | 5.0 |
| 50.0 | 0 | 100 | 0 | 5.0 |
| 50.1 | 100 | 0 | 0 | 5.0 |
| 60.0 | 100 | 0 | 0 | 0.1 |
| | | | | |

3 RESULTS

3.1 Standard Addition

HPLC analysis of solution A, the standard of 4-chloro-3-methylphenol, produced chromatogram A in figure 3.1 with a peak at 26.62 minutes. The 0.078 mg of 4-chloro-3-methylphenol that produced the peak in chromatogram A was from an injection volume of 100 μ L with a concentration of 0.780 mg/mL. Solution B, containing only MWF, produced a chromatogram with multiple peaks, but one was at 26.60 minutes. The

MWF sample solution contained 0.0751 g/ mL resulting in 7.51 mg of MWF being injected onto the column in chromatogram B. Solution C with standard addition resulted in chromatogram C with an increased height of the peak at 26.62 minutes. The increased peak height in solution C was the result of an added 0.082 mg or 0.82 mg/mL of 4-chloro-3-methlyphenol. Finally solution D with even more standard added resulted in an even higher peak area response at 26.62. See figure 3.1 The UV absorption peak area data for 4-chloro-3-methlyphenol from chromatograms a, b, and c were analyzed by least squares linear regression for their average relationship to the amount of 4-chloro-3-methlyphenol in the MWF. The linear equation found was $Y = 99734733 X + 8144831$. The X intercept at $Y = 0$ estimated the amount of 4-chloro-3-methlyphenol in the MWF at 0.082 mg. There was 7.512 mg of MWF injected. Therefore, there was 1.10% 4-chloro-3-methlyphenol (w/w) in the MWF. See figure 3.2

HPLC Chromatograms of 4-Chloro-3-methylphenol

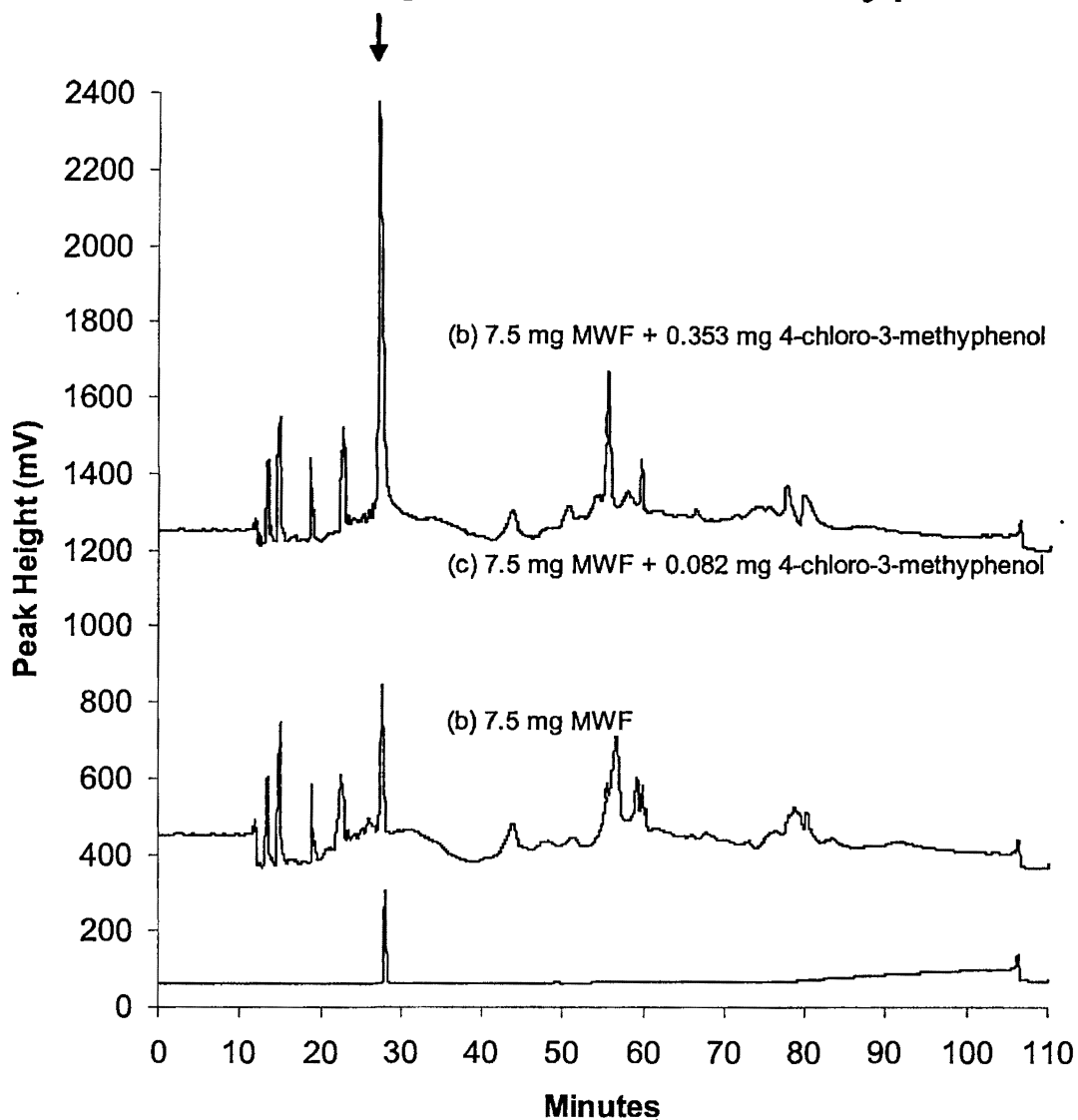


Figure 3.1

A) Chromatogram of 0.078 mg 4-chloro-3-methylphenol standard.

B) Chromatogram of 7.5 mg of MWF.

C) Chromatogram of 7.5 mg of MWF with 0.082 mg of 4-chloro-3-methylphenol added.

D) Chromatogram of 7.5 mg of MWF with 0.350 mg of 4-chloro-3-methylphenol added.

The retention time for 4-chloro-3-methylphenol was identified at 26.62 minutes from standard.

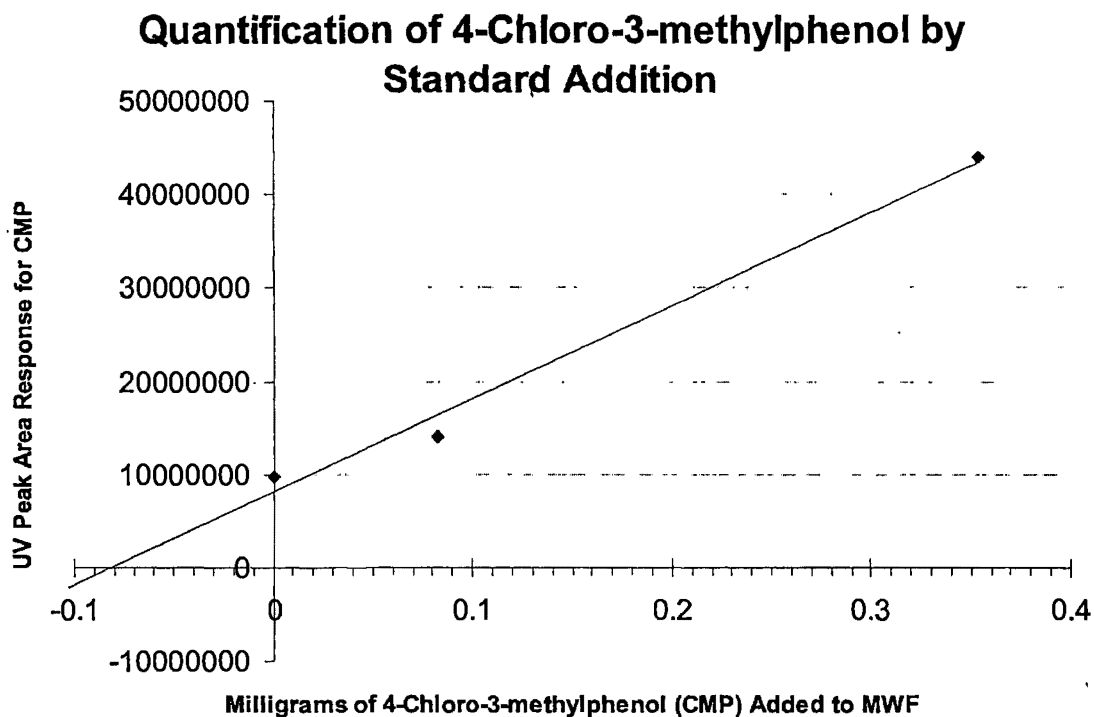


Figure 3.2 shows the quantification of 4-chloro-3-methylphenol by the standard addition technique. Least squares linear regression analysis of the 3 response data points against the amount of 4-chloro-3-methylphenol added resulted in a linear equation for peak area as a function of amount of analyte. The point where this line intercepts the X axis determined the amount of biocide in the MWF at 0.082 mg of 4-chloro-3-methylphenol in 7.512 mg of MWF or 1.10 %.

3.2 Fractionation

Figure 3.2 shows an overlay of the 17 chromatograms produced during the fractionation process. The two vertical lines that overlay the chromatograms intersecting the X axis at 27.0 and 29.0 represent the collection window for the 4-chloro-3-methylphenol peak. The mobile phase eluent containing 4-chloro-3-methylphenol was collected from each of the 17 runs between approximately 27.0 and 29.0 minutes and pooled. The mobile phase before and after the window was also collected and pooled. The two pooled fractions were concentrated using a nitrogen evaporator at 60 °C to a final volume of 1275 μ L back to the MWF samples original density of 1.0 g/mL.

Overlay of HPLC Fractionation Chromatograms

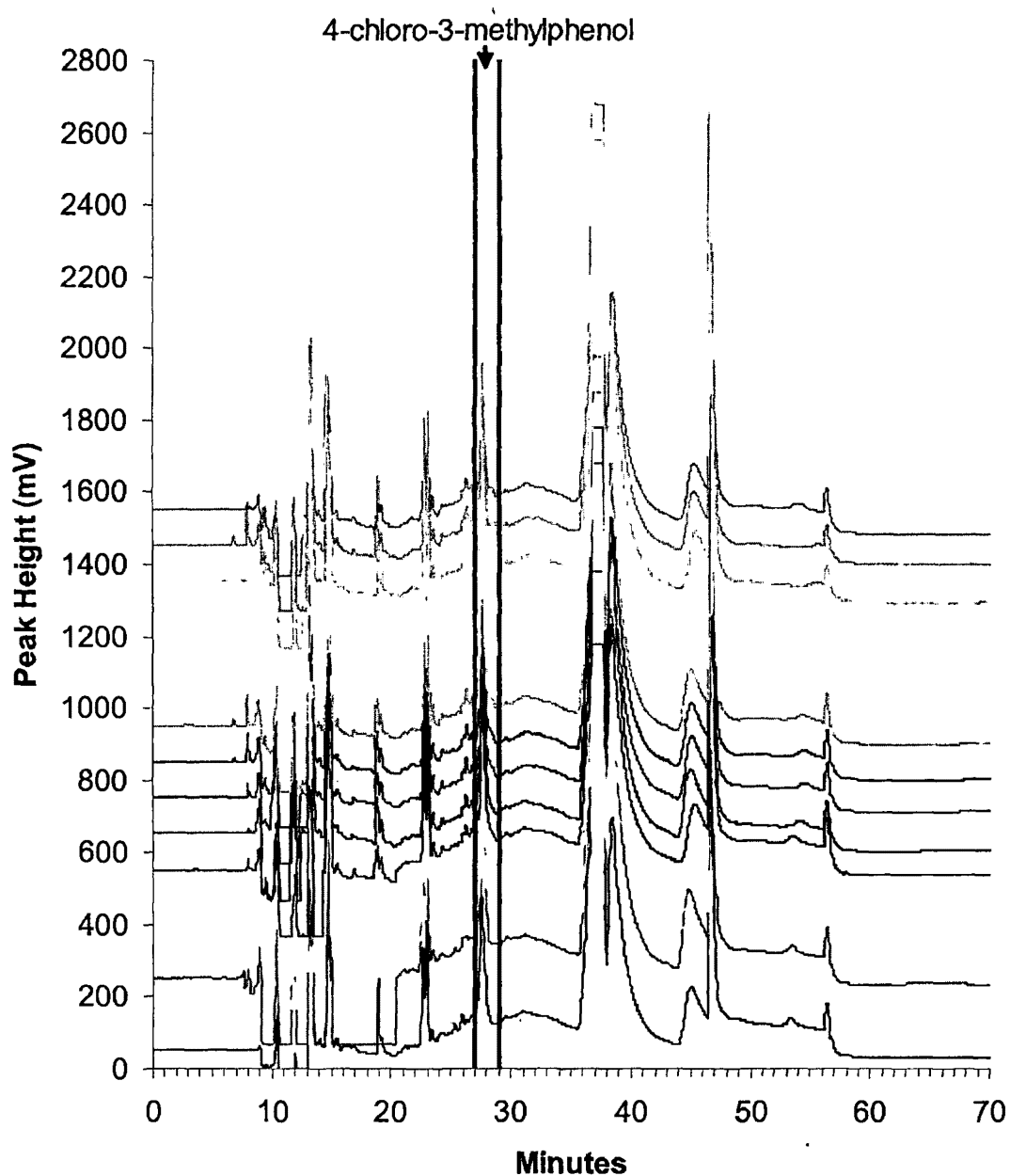


Figure 3.2

This is an overlay of the 17 runs. The runs remain in order with the first one being the bottom chromatogram. With exception to the first chromatogram, the peak heights for each run were offset by 100 units over the previous run. The two vertical lines represent the fractionation window of approximately 27.0 to 29.0 minutes.

3.3 COC-GC-MS Analysis of the Fractionations

The original MWF solution, the isolated 4-chloro-3-methylphenol fraction, and the MWF fraction minus the 4-chloro-3-methylphenol fraction were analyzed by COC-GC-MS. The total ion current (TIC) chromatogram resulting from COC-GC-MS analysis of the 4-chloro-3-methylphenol fraction showed only one peak after the solvent peak, Figure 3.3a. The identity of 4-chloro-3-methylphenol in that peak was confirmed by mass spectrometry through the high correlation of its mass spectrum with that of the library spectrum, Figure 3.3b. The parent ion at m/z of 142 was used in the future for more specific selected ion monitoring (SIM) analysis. An overlay of the three TIC chromatograms compared the 4-chloro-3-methylphenol content in each fraction Figure 3.3c. The 4-chloro-3-methylphenol fraction chromatogram showed its one peak, while the original MWF sample and MWF-4-chloro-3-methylphenol fractions showed similar matrix profiles but the matrix response masked whether or not it contained 4-chloro-3-methylphenol. Incidentally, a peak at about 18 minutes in the original MWF suffered loss during reconstitution. This peak was a pine scent component, terpineol, but otherwise the MWF fraction minus the 4-chloro-3-methylphenol appeared to contain all the other detectable components in the original matrix. An overlay of the more specific SIM chromatograms extracted from the TIC chromatogram, showed that only the original MWF and the 4-chloro-3-methylphenol fraction contained the biocide and that there was no trace of 4-chloro-3-methylphenol in the MWF fraction minus the 4-chloro-3-methylphenol fraction, i.e. it was completely resolved from the original MWF.

GC-MS Chromatogram of 4-chloro-3-methylphenol Fraction

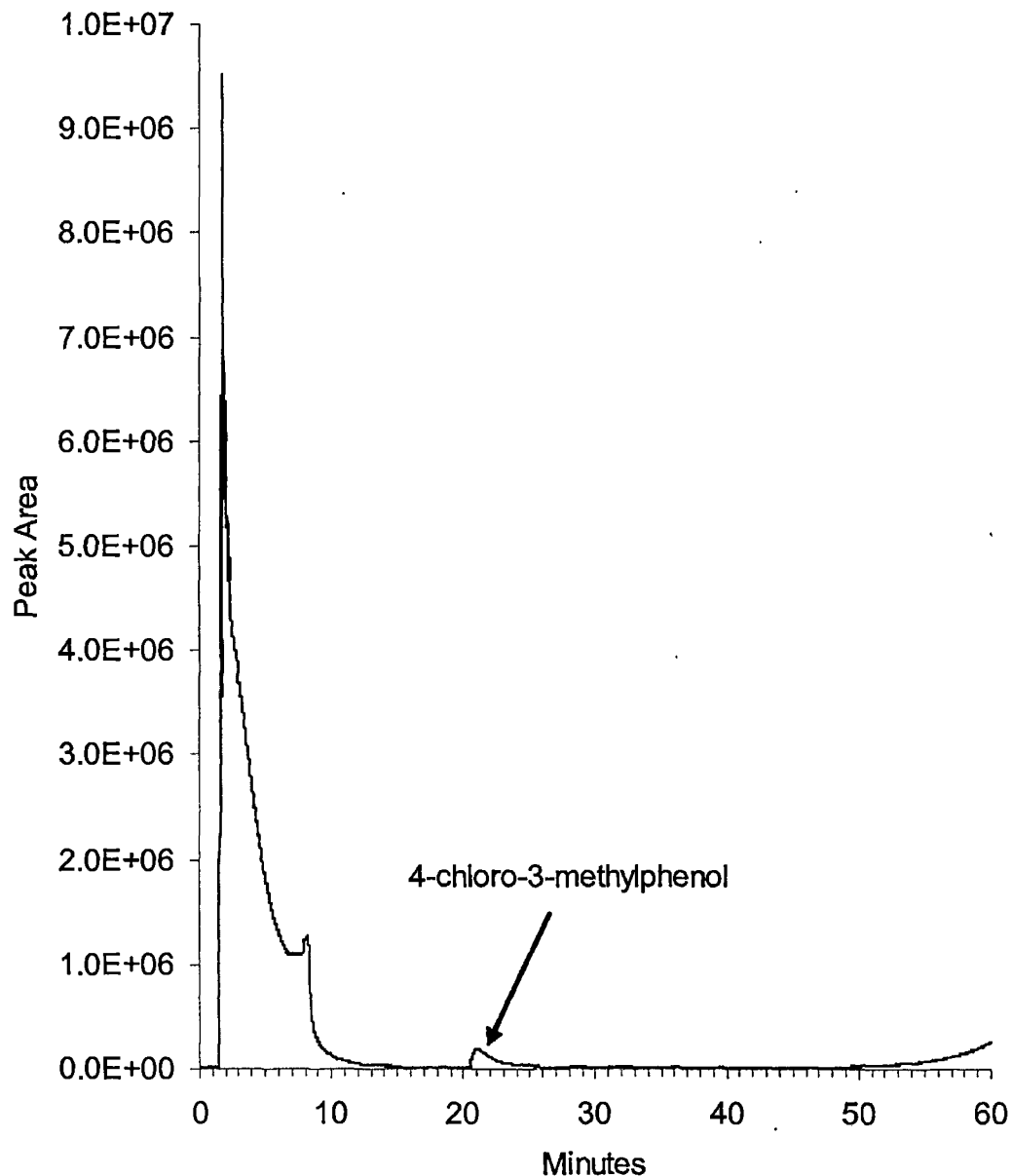


Figure 3.3a:

The suspect allergen was found in GC-MS chromatogram, the peak located at approximately 21 minutes.

**Mass Spectrum of 4-Chloro-3-methylphenol Fraction
Compared to Library Spectrum**

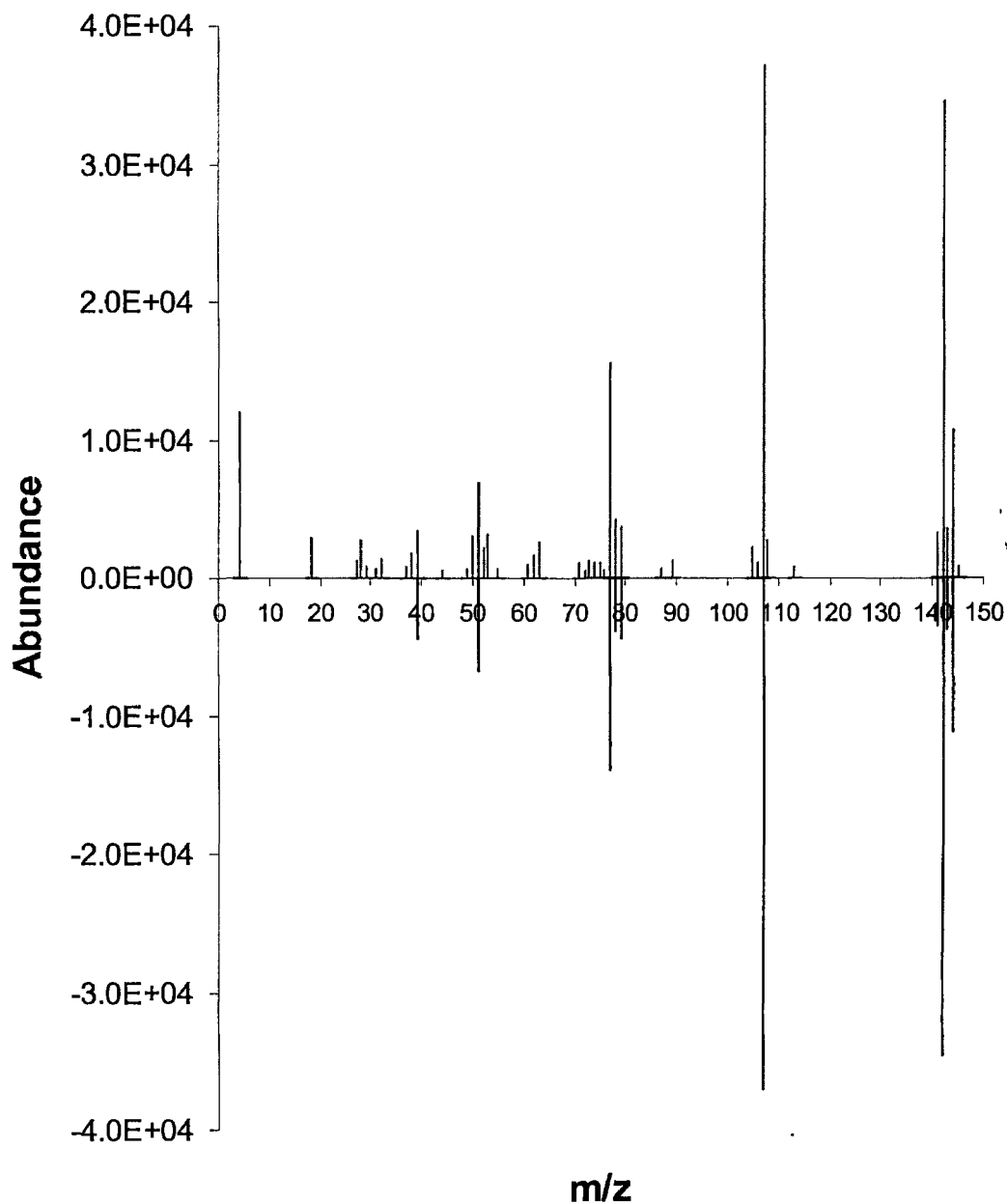


Figure 3.3b:

The mass spectrum of the 4-chloro-3-methylphenol peak at 21 minutes (top) is shown. The reference mass spectrum of 4-chloro-3-methylphenol (bottom) is from National Institute of Advanced Industrial Science and Technology (AIST)

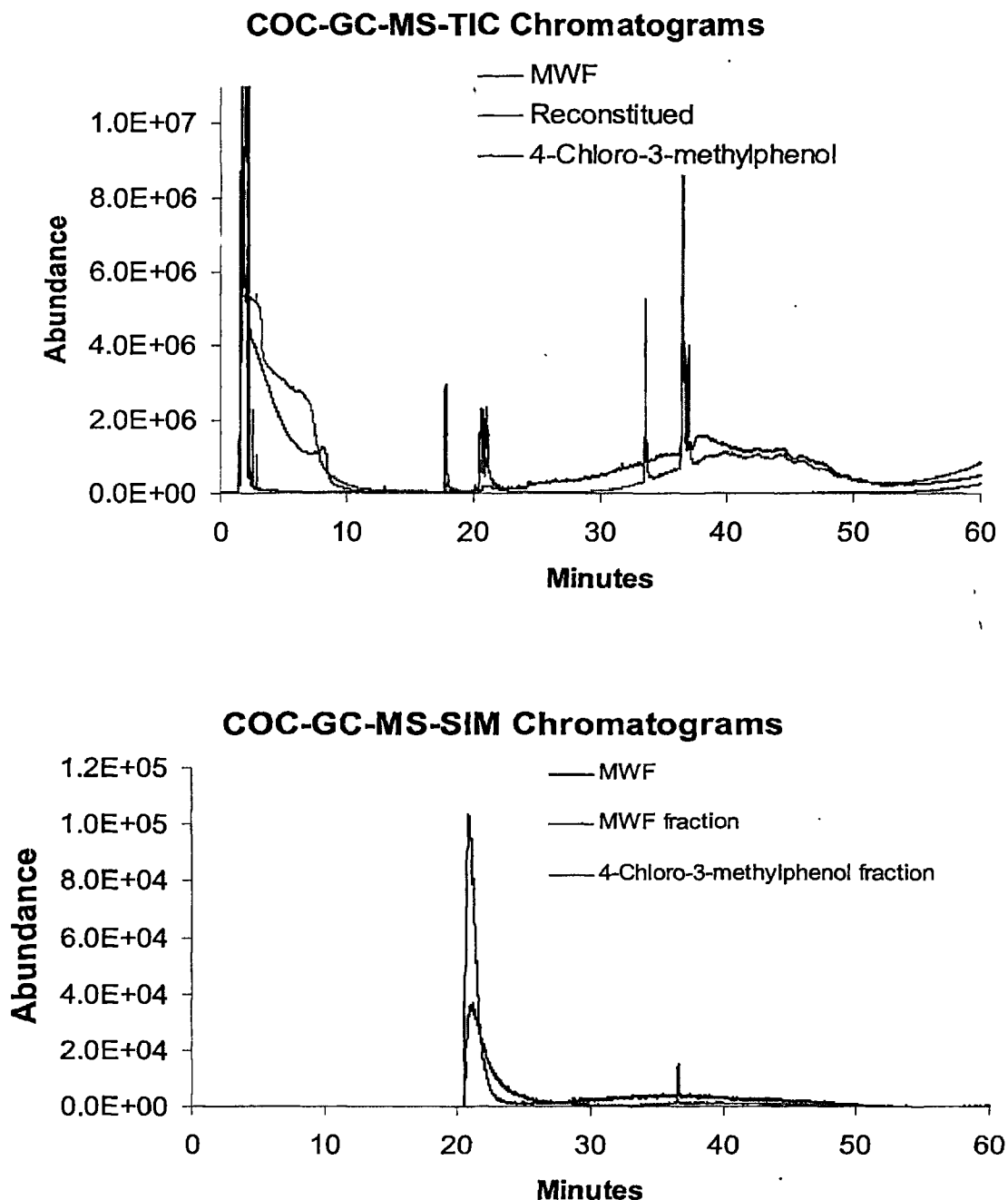


Figure 3.3c:

The chromatograms of the original MWF sample, the MWF fraction, and the 4-chloro-3-methylphenol (chlorocresol) fraction are shown. The upper chromatograms are of the sample using the cool-on-column injection gas chromatography with conventional electron impact mass spectrometry detection monitoring total ion current (COC-GC-MS-TIC). The lower chromatograms are same chromatograms, but with selective ion monitoring for the chlorocresol 142 m/z ion (COC-GC-MS-SIM).

4 CONCLUSION

An HPLC-UV method was developed that could separate, measure, and isolate the biocide, 4-chloro-3-methylphenol, in the MWF. The standard addition technique confirmed and quantified the content of 4-chloro-3-methylphenol in the MWF, and identified its retention time. COC-GC-MS analysis confirmed the presence of 4-chloro-3-methylphenol in the MWF, its absence in the MWF fraction, and its presence in the 4-chloro-3-methylphenol fraction. Multiple injections maximized the loading capacity of the semi-preparative column, 75 mg. It took 17 HPLC runs to fractionate 1275 mg of MWF for LLNA into two separated fractions. The semi-preparative HPLC method resulted in reproducible chromatograms for the seventeen fractionations of MWF. One run was sufficient to fractionate a metalworking fluid for COC-GC-MS analysis of the component. The two recovered fractions, the 4-chloro-3-methylphenol and the remaining MWF, were reconstituted by nitrogen evaporation to a volume of 1275 μL , i.e. to its original concentration. This method can be used to fractionate MWFs for COC-GC-MS analysis that requires trace enrichment and LLNA that requires 1000 mg of sample.

5 DISCUSSION

Previous studies have shown that this MWF was a potential cause of contact dermatitis and analysis of the fluid showed it to contain 4-chloro-3-methylphenol, a suspected allergen. Using HPLC, the 4-chloro-3-methylphenol was separated, identified, confirmed, quantified, and isolated from the MWF. This current analytical

method differs from IH methods where gravimetric analysis was used to determine the total mass concentrations of the MWF in air for assessment of worker exposure. In order to confirm that 4-chloro-3-methylphenol is an allergen, the LLNA needs to be performed on the separated fractions. If the fraction containing the suspect allergen produces the same allergenic response and the reconstituted MWF does not, then the allergen was identified.

6 SUGGESTIONS

To decrease the time it takes to process the desired amount of sample for LLNA, a larger HPLC preparative column and system can be used. The advantage of this system is that it used a standard HPLC analytical pump system, but the disadvantage was analytical pump systems are not designed for large preparative column flows and so the semi-preparative column was a compromise that required multiple runs. There are preparative HPLC systems commercially available that could have processed a gram in a single chromatogram. To decrease evaporation time, more than one evaporator can be utilized. Furthermore, this HPLC method appears to be versatile and applicable to many complex mixtures, and methodology may be used to identify and remove allergens from other complex mixtures. In the medical field specifically, HPLC in combination with the LLNA approach could be used to help identify and remove allergens that cause allergic contact dermatitis (ACD) in individuals.

REFERENCES

- [1] National Institute for Occupational Safety and Health. (2006).
<http://www.cdc.gov/niosh/topics/skin/skinresearch>
<http://www.cdc.gov/niosh/topics/metalworking>
- [2] Cheng, C., Phipps, D., & Alkhaddar, R. M. (2005). "Treatment of Spent Metalworking Fluids," *Water Research*, 39(17), 4051-63.
- [3] National Institute for Occupational Safety and Health. (1998). *NIOSH criteria for a recommended standard: Occupational exposure to metal working fluids* (DHHS (NIOSH) Publication No. 98-102). Cincinnati, OH: U.S. Department of Health and Human Services, Public Health Service, Centers for Disease Control and Prevention, National Institute for Occupational Safety and Health.
- [4] Sheehan, M. J. (1999). Final Report of the OSHA Metalworking Fluids Standards Advisory Committee. Occupational Safety and Health Administration (OSHA).
- [5] Lim, C. H., Yu, I. J., Kim, H. Y., Lee, S. B., Kang, M. G., Marshak, D. R., & Moon, C. K. (2005). Inflammatory and Immunological Responses to Subchronic Exposure to Endotoxin-Contaminated Metalworking Fluid Aerosols in F344 Rats. *Environmental Toxicology*, 20, 212-218.
- [6] Glaser, R., Pretty, J., Brown, K., Arnold, J., Lunsford, R. A., & Park, S.H. (2003). "Analytical method research for metalworking fluids at the National Institute for Occupational Safety and Health (NIOSH)," *Gefahrstoffe – Reinhaltung der Luft*, 62, 237-240.
- [7] National Institute for Occupational Safety and Health. (2003). *NIOSH Manual of Analytical Methods (NMAM)*, 4th ed. (DHHS (NIOSH) Publication No. 94-113). Washington, DC: U.S. Department of Health and Human Services, Public Health Service, Centers for Disease Control and Prevention, National Institute for Occupational Safety and Health.
- [8] Verma, D. K., Shaw, D. S., Shaw, M. L., Julian, J. A., McCollin, S. A., & des Tomber, K. (2006), "An Evaluation of Analytical Methods, Air Sampling Techniques, and Airborne Occupational Exposure of Metalworking Fluids," *Journal of Occupational and Environmental Hygiene*, 3, 53-66.
- [9] Raynor PC, Kim SW, & Bhattacharya M. (2005), "Mist generation from metalworking fluids formulated using vegetable oils." *Ann Occup Hyg.*, 49(4):283-93
- [10] Lee, B. L., Ong, H. Y., Koh, D., and Ong, C. N. (1994), "High-performance liquid chromatographic method for determination of dehydroabietic and abietic acids, the skin sensitizers in binding adhesive," *Journal of Chromatography A*, 685, 263-269.

- [11] Wahl, R., Meineke, D., & Maasch, H. J. (1987). Molecular Weight Determination of Allergen Extracts and Isolation of Allergenic Molecules by High-Performance Liquid Chromatography. *Journal of Chromatography A*, 397, 307-312.
- [12] Andersen K. E., Carlsen L., Egsgaard H, & Larsen E. (1985), "Contact sensitivity and bioavailability of chlorocresol," *Contact Dermatitis* 13 (4), 246–251.
- [13] Final Report on the Safety Assessment of p-Chloro-*m*-cresol, *International Journal of Toxicology*, 16(3), 235-268.

Society of Tribologists & Lubrication Engineers

**62nd Annual Meeting of the
Society of Tribologists and
Lubrication Engineers
2007**

**May 6-10, 2007
Philadelphia, Pennsylvania, USA**

Volume 1 of 2

Printed from e-media with permission by:

Curran Associates, Inc.
57 Morehouse Lane
Red Hook, NY 12571
www.proceedings.com

ISBN: 978-1-60423-326-1

Some format issues inherent in the e-media version may also appear in this print version.