
Job Aid for Clinicians

How to send information to the health department about a patient under investigation (PUI) for AFM



Contact your health department when you identify a PUI for AFM.



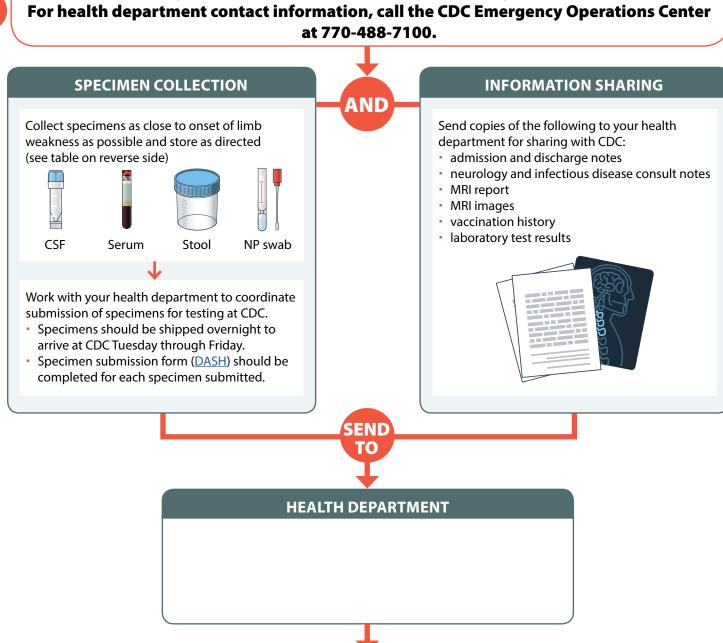



Specimens to collect and send to CDC for testing for AFM PUIs

SAMPLE	AMOUNT	TUBE TYPE	PROCESSING	STORAGE	SHIPPING
CSF	1mL (collect at same time or within 24hrs of serum)	Cryovial	Spun and CSF removed to cryovial	Freeze at -70°C	Ship on dry ice
Serum	≥0.4mL (collect at same time or within 24 hours of CSF)	Tiger/red top	Spun and serum removed to tiger/red top.	Freeze at -70°C	Ship on dry ice
Stool	≥1 gram (2 samples collected 24hrs apart)	Sterile container	n/a	Freeze at -20°C *	Ship on dry ice. Rectal swabs should not be sent in place of stool.
Respiratory (NP)/ Oropharyngeal (OP) swab	1ml (minimum amount)	n/a	Store in viral transport medium	Freeze at -20°C *	Ship on dry ice

^{*}All specimens may be stored at -70 $^{\circ}$ C for ease of shipping.

Coordinate with your health department to send information about PUIs and ship specimens to CDC. For health department contact information, call the CDC Emergency Operations Center at 770-488-7100.

